


Trelstar®

Instructions for Use

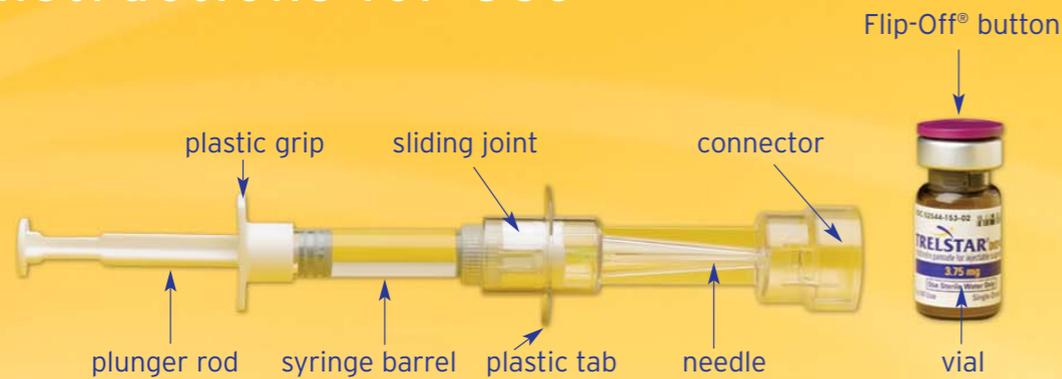

Clip'n'Ject® Preparation

Wash your hands with soap and hot water and put on gloves immediately prior to preparing the injection. Place the package containing the Clip'n'Ject system and the Trelstar® vial on a clean, flat surface that is covered with a sterile pad or cloth. Peel the Tyvek® cover away from the blister package, and place the vial, connector, alcohol swab, and plunger rod on the prepared surface.



1. Be sure to begin by removing the Flip-Off® button from the top of the vial, revealing the rubber stopper.



2. Disinfect the rubber portion of the vial cap with the alcohol swab. Discard the alcohol swab and let the alcohol dry.

3. Proceed to Clip'n'Ject Activation.

Clip'n'Ject® Activation




1. Holding the vial upright and flat on the table surface with one hand, place the plastic connector directly over the top of the Trelstar® vial with the other hand. Press the connector down firmly on the vial top. This will ensure proper positioning of the vial.




2. Still holding the vial with one hand, press the syringe barrel downward as far as it will go in the connector. This results in insertion of the needle into the rubber stopper in the vial top to the predetermined depth.



3. Check to make sure that the needle is inserted into the vial. Now, screw the plunger rod into the end of the plastic grip on the syringe barrel.

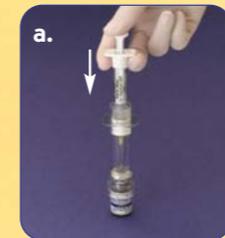

4. **a.** With the vial still on the flat surface, place your thumb on the plunger rod and depress the plunger rod to inject the sterile water diluent into the vial. **b.** With your thumb on the plunger rod, place two fingers under the plastic tab on the connector to keep the assembly together. Gently rotate the system so that the diluent rinses the vial sides to ensure complete mixing of Trelstar® and the sterile water diluent. The solution will now have a milky appearance. In order to avoid separation of the solution, proceed to the next steps without delay.

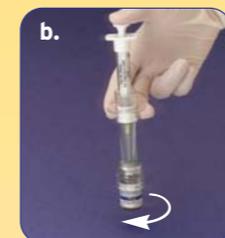

5. Hold the Clip'n'Ject system in a vertical position with the connector at 12 o'clock and the syringe plunger rod at 6 o'clock. Double check to make sure that the syringe is still as far forward as possible in the connector with the needle situated in the vial. Grasp the Clip'n'Ject system firmly by the syringe barrel and pull back the plunger rod to draw the reconstituted Trelstar® into the syringe.



6. Immediately before injecting Trelstar®, remove the filled syringe from the connector by holding the syringe by the barrel and pressing your thumbs against the plastic tabs of the connector and pulling the syringe section from the connector. Trelstar® is now ready for administration.

Clip'n'Ject® Disposal



After administering Trelstar®, dispose of the Clip'n'Ject system as follows:

- a.** Place Clip'n'Ject with attached vial in standing upright position on a flat surface.
- b.** Using one hand, replace the syringe into the Clip'n'Ject connector.
- c.** Dispose of syringe and attached Clip'n'Ject connector with vial into a suitable sharps container.

TRELSTAR[®] LA
11.25 mg
(triptorelin pamoate for injectable suspension)

TRELSTAR[®] DEPOT
3.75 mg
(triptorelin pamoate for injectable suspension)

Clip'n'Ject and Flip-Off are registered trademarks of West Pharmaceutical Services, Inc. Tyvek is a registered trademark of E.I. du Pont de Nemours and Company. © 2005, Watson Pharma, Inc., Morristown, NJ 07962. All rights reserved. 03493 9/05

WATSON Pharma, Inc.
A Subsidiary of Watson Pharmaceuticals, Inc.